

Thompson Retractor

Mini Retractor Lite LED


Replace your headlight

Improve your exposure beyond retraction by adding cool, bright light to your cervical spine procedure with Thompson's low profile Mini Retractor Lite LED.

The Mini Retractor Lite LED is indicated for use in cervical procedures only.




Utilize the bright, reusable Mini Retractor Lite LED to add illumination to your cervical procedures.

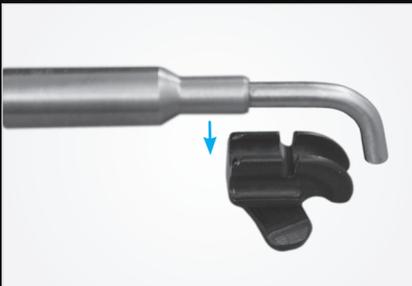

Step 1

Press distal end of Mini Retractor Lite LED (#40002MRL) into Lite Clip (#40026MRL).

Apply pressure to the rigid metal portion of the light only to avoid damaging cable.

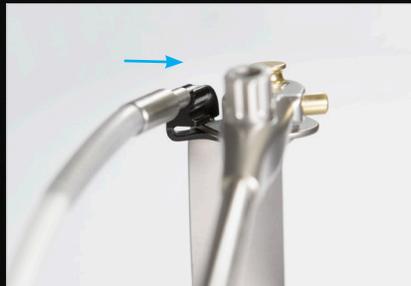

Step 2

Slide onto retractor blade (light is intended to be used with Thompson cervical blades, see Mini Retractor Lite LED IFU for exceptions).

Adjust clip so that the light is not obstructed and points down blade face.

Plug cable into LED Light Source (#LLS-2000-2) and turn source on.



Tip

Avoid moving Mini Retractor Lite cable and reduce cable tension to prevent clip from moving on the blade.



NOTICE:

S-Lock may be utilized to rotate blade without affecting clip position in most cases. If clip position is affected, move clip and light to a different blade.



NOTICE:

See Mini Retractor Lite LED IFU for cleaning, sterilization, and care instructions.



MINI RETRACTOR LITE LED KIT * #40001MRL

REF	QTY	ITEM DESCRIPTION	PART #
A	1	Mini Retractor Lite LED	40002MRL
B	2	Mini Retractor Lite Clip (Radiopaque)	40026MRL

REF	QTY	ITEM DESCRIPTION	PART #
	1	Instrument Case (11" x 7" x 3 1/2")	50000T
C	1	LED Light Source	LLS-2000-2

* Cervical use only



NOTICE:

LED Light Source (#LLS-2000-2) required for use.

Thompson 
 SURGICAL INSTRUMENTS
 10341 East Cherry Bend Road
 Traverse City, Michigan 49684
 phone: 231.922.0177
 fax: 231.922.0174
thompsonretractor.com

* Free trial valid for U.S. customers only. International customers please call for more information.

© 2024 Thompson Surgical Instruments, Inc. Traverse City, Michigan.

® S-Lock®, PLA®, and the "T Circle" logomark  are Registered Trademarks of Thompson Surgical Instruments, Inc.

Patents: US4971038, US5025780, US5888197, US5897087, US5902233, US5984865, US6033363, US6416465, US6511423, US7338442, US7749163, US8257255, US8360971, US8617064

Other patents pending.

The Mini Retractor Lite LED (40002MRL) and LED Light Source (LLS-2000-2) are manufactured by Sunoptic Technologies® and distributed by Thompson Surgical Instruments.