


"Thompson's retractors are held more stably and have more versatility than an assistant."

- Jeffrey P. Beckenbaugh, DO

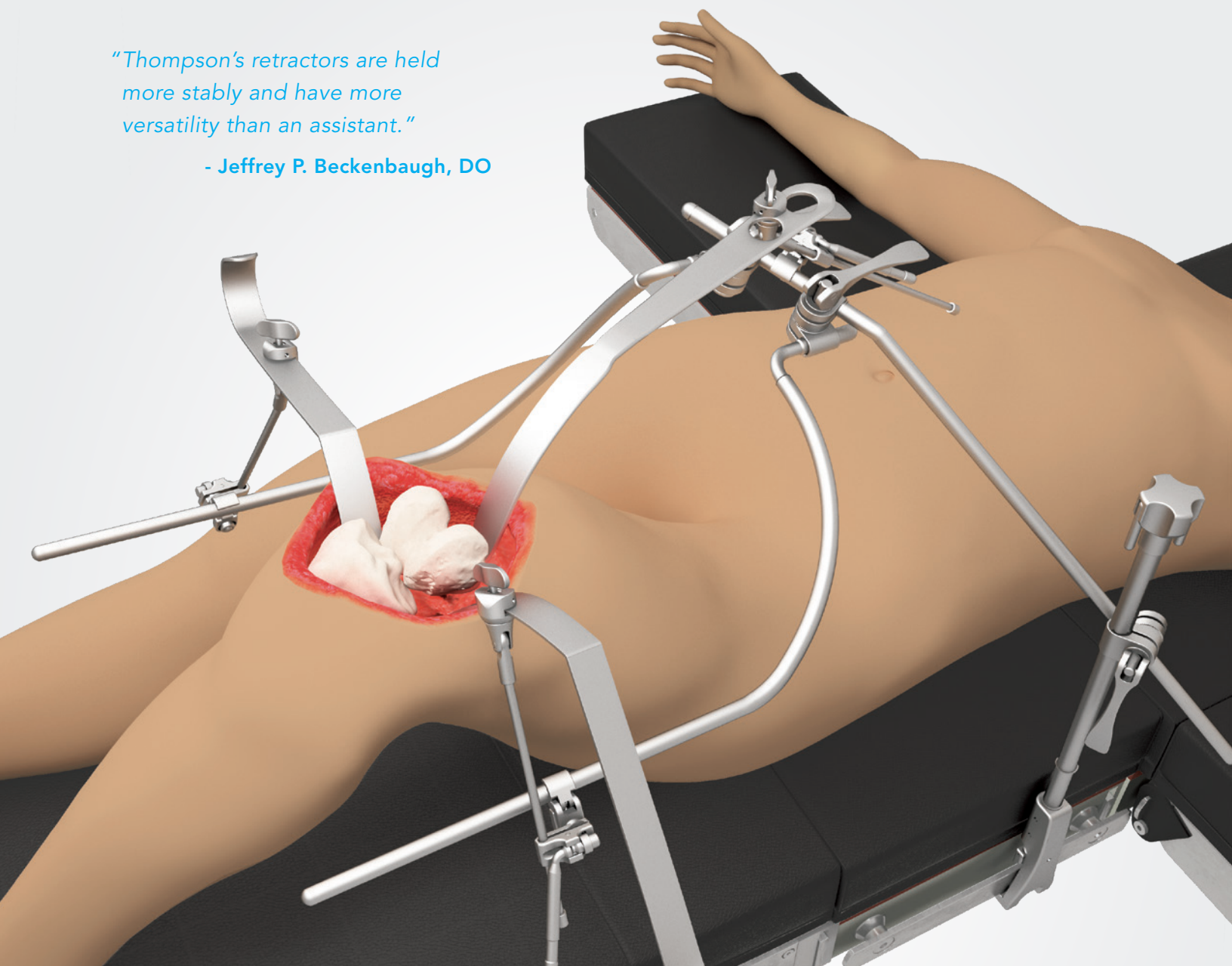

Uncompromised Exposure™ for **Knee** with the **Thompson Retractor Total Joint System**

Table-Mounted Total Joint Retraction


Firm Hold with Fewer Hands Required

Thompson's table-mounted Total Joint Retractor System is ideal for orthopedic surgeons performing hip, knee, and shoulder procedures. This specialized system utilizes rigid handles or flexible tethers to attach retractor blades to a stable, table mounted frame which eliminates the need for extra staff to hold retractors. A sterile mount, low profile design, and streamlined set up saves time and space.

Thompson Retractor Exclusive Features: Thompson is dedicated to offering the most forward thinking, proven, and high quality features. These exclusive, often patented, features define Thompson and allow for Uncompromised Exposure.

"The Thompson Retractor makes my assistant's job a breeze and is invaluable in the difficult or stiff knee."

- Jeffrey P. Beckenbaugh, DO

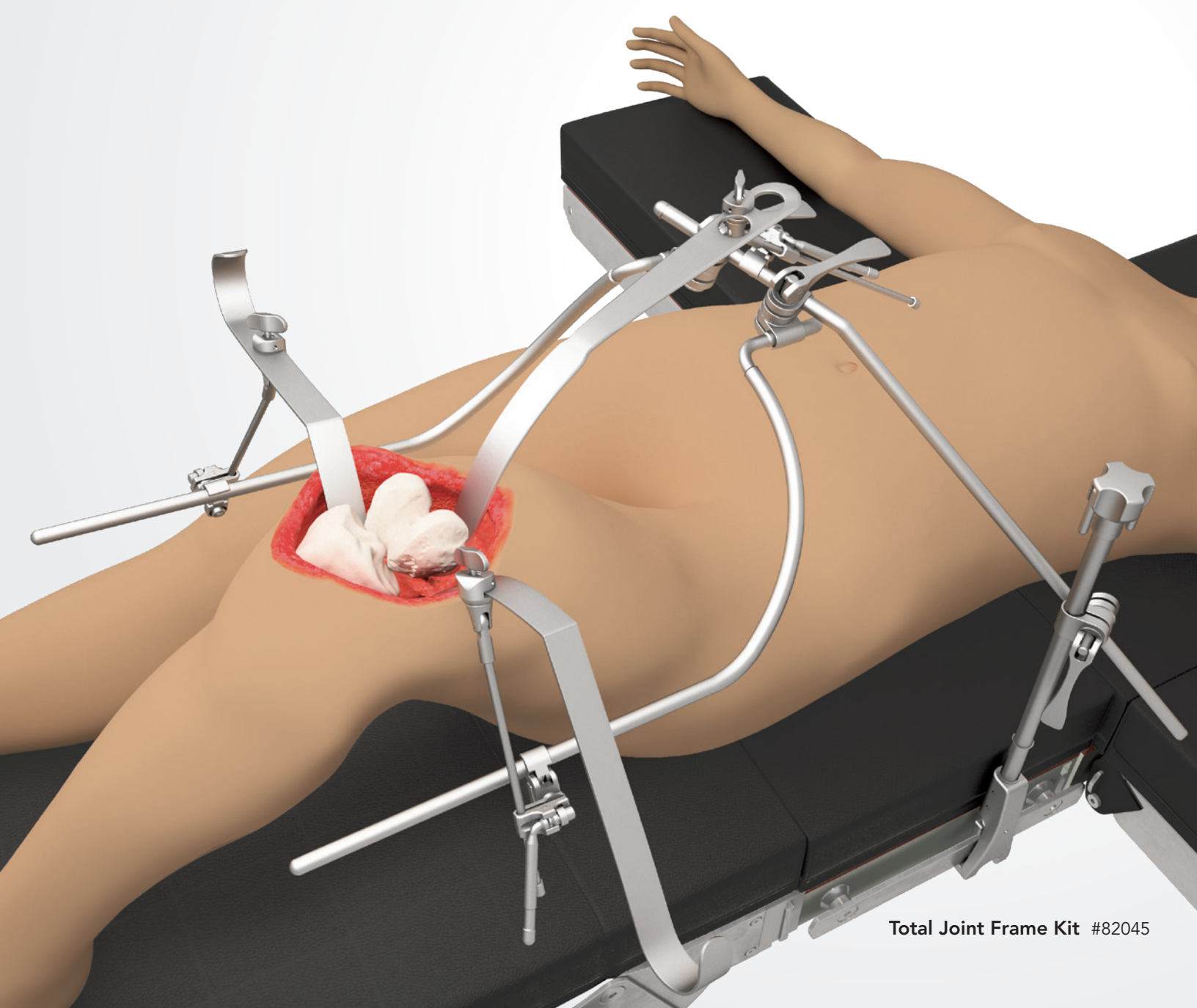

Total Joint Frame Kit #82045



Compatibility

Thompson's Total Joint Retractor System is compatible with most orthopedic retractor blades.



Anatomically Designed Frame

Low profile, independent frame arms contour to patient anatomy.



Sterile Mount Rail Clamp

Thompson's table mounted frame eliminates retractor shifting. Top control allows attachment to table, as well as adjustments, to be achieved within the sterile field.



Hands-Free Exposure

Photo courtesy Jeffrey P. Beckenbaugh, DO

Uncompromised Exposure™



Tethers

Flexible tethers attach quickly and securely while remaining flexible to allow slight movements when needed.

#42105 (5-pack)



Universal Blade Adapter Handle

Adapter Handles attach quickly along flat blade surface, providing a secure, rigid hold.

#42126BP



Orthopedic Clip-On Hook Handle

Fast, easy option for blade attachment and adjustments. Compatible with most hip blades containing holes.

#42126OHP

For a Free Trial Call Today*

1.800.227.7543

Rev A
030922
mg tjksq



Uncompromised Exposure™

VISIT US ONLINE: THOMPSONSURGICAL.COM

CONTACT

Address

Thompson Surgical Instruments, Inc.
10341 E Cherry Bend Rd.
Traverse City, MI 49684
USA

Phone + Fax

Toll Free 1-800-227-7543
Phone 231-922-0177
Fax 231-922-0174

Online

Email contact@thompsonsurgical.com
Website www.thompsonsurgical.com

Social



Thompson
SURGICAL INSTRUMENTS

**Made In
USA**

* Free trial valid for U.S. customers only. Customers outside the U.S.A., please call +1-231-922-0177 for availability.

© 2022 Thompson Surgical Instruments, Inc. Printed in U.S.A.

© S-Lock®, PLA®, and the "T Circle" logomark are Registered Trademarks of Thompson Surgical Instruments, Inc.

Patents: trpat.com

Note: As we continually strive to provide the best products possible, some of the images in this document may appear slightly different from the product received.