



Carotid

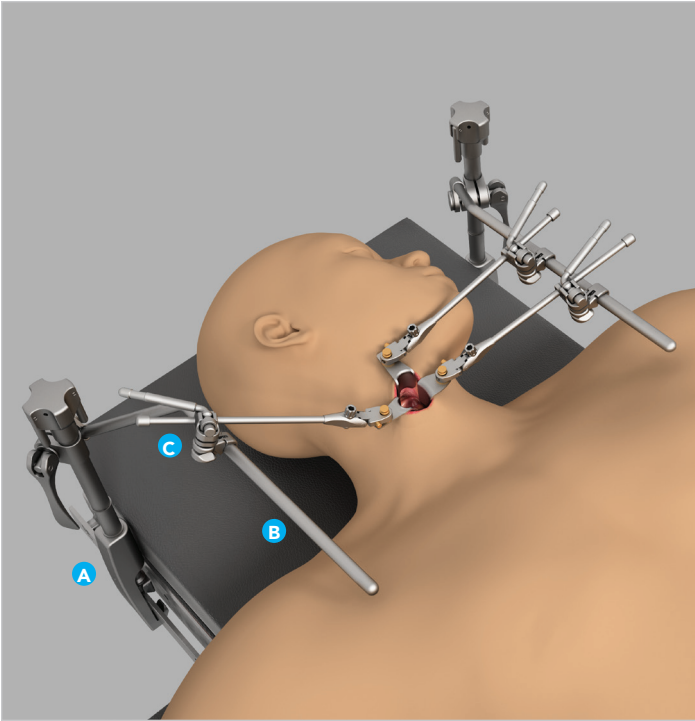

Thyroid

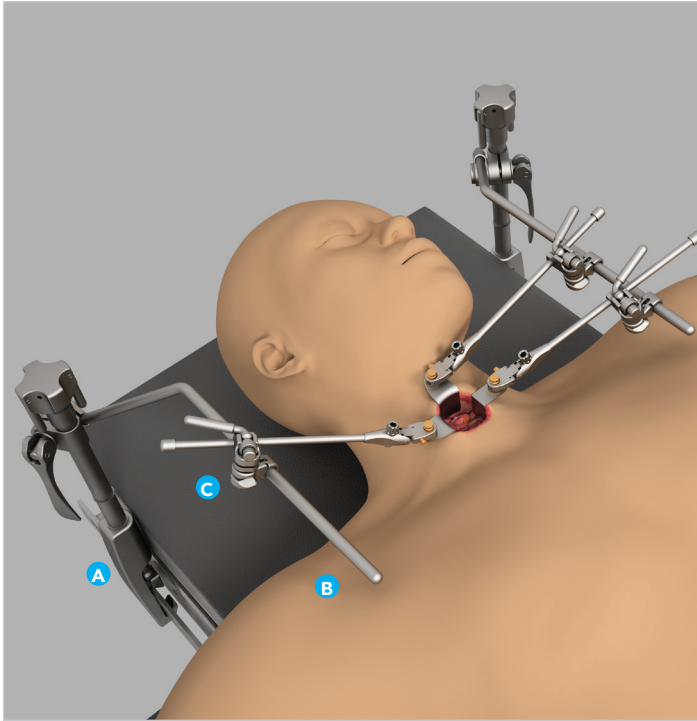

SMALL INCISION SYSTEM #SL82037

Non S-Lock system also available, may require a longer lead time.

REF	IMAGE	PART #	ITEM DESCRIPTION	QTY	APPLICATION
Tray 1					
A		43902ACS	Elite III Rail Clamp with 1 Cam Joint 11 "	2	Attaches to the rail of the OR table and holds the angled arm in position.
B		44117N	17" Angled Arm (10" x 7") @ 45°	2	
		44120	20" Angled Arm (10" x 10") @ 45°	2	Attaches to the rail clamp and holds the retractor handles in place around the incision.
		44012H	12" Extension Arm with Spring Hinge Joint	1	
C		SL42126WGP	Cam II Clip-on Quick Angle 8"	5	Attaches retractor blades to frame and angles retractor blades as needed, 45° up or down.
		60020	T-Handle	2	May be used on Quick Angle Handle to angle blade 45° up or down.

SMALL INCISION SYSTEM #SL82009 (Continued)

REF	IMAGE	PART #	ITEM DESCRIPTION	QTY	APPLICATION
		SL46125	Richardson 19mm x 25mm ($\frac{3}{4}$ " x 1")	2	Muscle retraction for exposure with small or mini-incision.
		SL46125A	Richardson 19mm x 51mm ($\frac{3}{4}$ " x 2")	1	
		SL46126	Richardson 25mm x 32mm (1" x 1 $\frac{1}{4}$ ")	1	
		SL46149	Richardson 44mm x 51mm (1 $\frac{3}{4}$ " x 2")	2	Abdominal wall retraction.
		SL46150	Richardson 51mm x 102mm (2" x 4")	1	
		SL46115	Malleable 13mm x 64mm ($\frac{1}{2}$ " x 2 $\frac{1}{2}$ ")	2	Soft tissue retraction.
		SL46211	Volkman Blunt 6 Prong 25mm x 13mm (1" x $\frac{1}{2}$ ")	1	Skin and shallow muscle retraction.
		SL46212	Volkman Sharp 4 Prong 16mm x 13mm ($\frac{5}{8}$ " x $\frac{1}{2}$ ")	1	
		SL46213	Volkman Blunt 4 Prong 16mm x 13mm ($\frac{5}{8}$ " x $\frac{1}{2}$ ")	1	
		SL46214	Murphy Sharp 4 Prong 25mm x 13mm (1" x $\frac{1}{2}$ ")	1	
		SL46215	Murphy Blunt 4 Prong 25mm x 13mm (1" x $\frac{1}{2}$ ")	1	
		SL46217	Murphy Blunt 3 Prong 19mm x 13mm ($\frac{3}{4}$ " x $\frac{1}{2}$ ")	1	
		SL45157A	Concave Lateral 25mm x 40mm (1" x 1 $\frac{1}{2}$ ")	1	Medial and lateral retraction for head and neck procedures.
		SL45159A	Concave Lateral 25mm x 70mm (1" x 2 $\frac{3}{4}$ ")	2	
		SL45232A	Concave Medial 25mm x 50mm (1" x 2")	2	Medial retraction for head and neck procedures.
		SL45153A	Concave Medial 32mm x 70mm (1 $\frac{1}{4}$ " x 2 $\frac{3}{4}$ ")	2	
		50000G	Instrument Case (23" x 11" x 3 $\frac{1}{2}$ ")	1	Holds instruments for sterilization, storage, and transportation.